

PATIENT CLINICAL INFORMATION

Patient's Name

Practitioner

41998318 NAME: CONSERVA,MAX FABRAZI
DOB: 02/02/1981 M COPAY:
FSC: 743 PLAN: ANTHEM BX/BS OF MI 210 2S
ATTG MD: UNKNOWN, SCHED 99996
REF MD: NO, REFERRING 99997
AUTH#
VISIT: 19603373

STOR: OAF LOC: PS9
DOS: 08/22/2011

Added a 1/2" thick
XPE pad to reduce

external rotation of femur in KAFO.
Pt to try + return later to feedback.
Gave him cash estimate.

~~Panel~~ CO

41998318 NAME: CONSERVA,MAX FABRAZI
DOB: 02/02/1981 M COPAY:
FSC: 743 PLAN: ANTHEM BX/BS OF MI 210 2S
ATTG MD: UNKNOWN, SCHED 99996
REF MD: NO, REFERRING 99997
AUTH#
VISIT: 19753620

STOR: OAF LOC: PS9
DOS: 09/23/2011

Moved pad more
posteriorly to ↓ rubbing
directly against femur + catching it

posteriorly. Also added aliplast to it +
gave pt extra gel cover, shearbar + velcro.
Ret in 1 week.

~~Panel~~ CO

PATIENT CLINICAL INFORMATION


Patient's Name

Practitioner

8/10/11 (cont) KAFO he feels pushes his prox tibia forward as he bends (max flex $\sim 70^\circ$) \rightarrow I think we can correct this \bar{c} Townsend free motion steel knee joint, $\sim 3/16" \times 5/8"$.

also complains of poor thigh internal rotation control as he flexes, and external rotation as his knee extends \rightarrow dig into thigh to hold femur, similar to hi fidelity modifiers.

I'll pad current thigh shell upon next visit. I will give pt \$ estimate for a new KAFO.

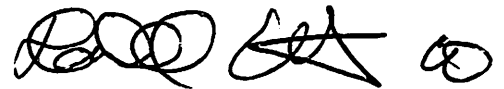

PATIENT CLINICAL INFORMATION

Patient's Name

Practitioner

41998318 NAME: CONSERVA, MAX FABRAZI
 DOB: 02/02/1981 M COPAY:
 FSC: 743 PLAN: ANTHEM BX/BS OF MI 210 2S
 ATTG MD: UNKNOWN, SCHED 99996
 REF MD: NO, REFERRING 99997
 AUTH# STOR: OAF LOC: PS9
 VISIT: 19493595 DOS: 07/28/2011

Pt is old MVA ^{in 1989,} extensive trauma to RLE.

Lost much of his lateral loose ^{+ ligaments} at R knee → severe g. valgus instability. Walt + Julie made a carbon fiber KAFO with multiple modifications that he's used for ~3 yrs. Needs ankle strap repaired. I did this. Basically very pleased ^{with this KAFO.}

Pt wants new KAFO, carbon thigh is solid ankle plastic. Wants more lateral thigh + calf control than current KAFO. Current KAFO has bilateral free ROM - knee joints, but thigh + calf shells overlap at knee for comfortable control at medial knee - just needs more of shells medially + needs plastic ^{in BK section}. Discuss \$ + cast later.

41998318 NAME: CONSERVA, MAX FABRAZI
 DOB: 02/02/1981 M COPAY:
 FSC: 743 PLAN: ANTHEM BX/BS OF MI 210 2S
 ATTG MD: UNKNOWN, SCHED 99996
 REF MD: NO, REFERRING 99997
 AUTH# STOR: OAF LOC: PS9
 VISIT: 19546565 DOS: 08/10/2011

~~Handwritten signature~~

Pt's R LE is 1 5/8" shorter than

L. Order ShearBan. R patella is rotated externally ~30° + ~1" too high. Current

PATIENT CLINICAL INFORMATION

Patient's Name

Practitioner

41998318 NAME: CONSERVA,MAX FABRAZI
DOB: 02/02/1981 M COPAY:
FSC: 743 PLAN: ANTHEM BX/BS OF MI 210 2S
ATTG MD: UNKNOWN,SCHED 99996
REF MD: NO,REFERRING 99997
AUTH# STOR: OAF LOC: PS9
VISIT: 16402214 DOS: 08/04/2009

41998318 NAME: CONSERVA,MAX FABRAZI
DOB: 02/02/1981 M COPAY:
FSC: 743 PLAN: ANTHEM BX/BS OF MI 210 2S
ATTG MD: UNKNOWN,SCHED 99996
REF MD: NO,REFERRING 99997
AUTH# STOR: OAF LOC: PS9
VISIT: 16480764 DOS: 08/21/2009

41998318 NAME: CONSERVA,MAX FABRAZI
DOB: 02/02/1981 M COPAY:
FSC: 743 PLAN: ANTHEM BX/BS OF MI 210 2S
ATTG MD: UNKNOWN,SCHED 99996
REF MD: NO,REFERRING 99997
AUTH# STOR: OAF LOC: PS9
VISIT: 16561022 DOS: 09/11/2009

Fig 8 ankle strap replaced. Pt is comfortable w/ today's
adjustment and will contact further per



41998318 NAME: CONSERVA,MAX FABRAZI
DOB: 02/02/1981 M COPAY:
FSC: 743 PLAN: ANTHEM BX/BS OF MI 210 2S
ATTG MD: UNKNOWN,SCHED 99996
REF MD: NO,REFERRING 99997
AUTH# STOR: OAF LOC: PS9
VISIT: 19796446 DOS: 09/28/2011

PATIENT CLINICAL INFORMATION

Patient's Name

Practitioner

41998318 NAME: CONSERVA,MAX FABRAZI
DOB: 02/02/1981 M COPAY: 15
FSC: 743 PLAN: BLUE CROSS/SHLD OF MI 2SU
ATTG MD: UNKNOWN,SCHED 99996
REF MD: NO,REFERRING 99997
AUTH# STOR: OAF LOC: PS9
VISIT: 15806326 DOS: 03/20/2009

3/20/09

pt in to determine if we could do
a new KAPO on his current COBRA insurance
which ends April 15. will find out and contact
him.

Julie J. C.O.

3/23/09

Spoke with max on the telephone. Unable to
do on current insurance. will wait a few months
and try on new insurance.

Julie J. C.O.

PATIENT CLINICAL INFORMATION

Patient's Name

Practitioner

41998318 NAME: CONSERVA, MAX FABRAZI
DOB: 02/02/1981 M COPAY: 15
FSC: 743 PLAN: BLUE CROSS/SHLD OF MI 2SU
ATTG MD: UNKNOWN, SCHED 99996
REF MD: DELGADO, ELIANA D MD 34842
AUTH# STOR: OAF LOC: P59
VISIT: 15395516 DOS: 12/12/2008

Pt enjoys new AFO piece. needs to be moved 1/4" distally. Did that. pt pleased

Julia

41998318 NAME: CONSERVA, MAX FABRAZI
DOB: 02/02/1981 M COPAY:
FSC: 743 PLAN: BLUE CROSS/SHLD OF MI 2SU
ATTG MD: UNKNOWN, SCHED 99996
REF MD: NOT AVAILABLE 99995
AUTH# STOR: LOC: HLD
VISIT: 15425297 DOS: 12/16/2008

lateral malleolus has pressure from carbon section that held. cut off extension. added rivets to reduce motion between AFO and calf piece. Will return after X-mas to discuss new KAFO.

Julia

41998318 NAME: CONSERVA, MAX FABRAZI
DOB: 02/02/1981 M COPAY: 15
FSC: 743 PLAN: BLUE CROSS/SHLD OF MI 2SU
ATTG MD: UNKNOWN, SCHED 99996
REF MD: NOT AVAILABLE 99995
AUTH# STOR: LOC: HLD
VISIT: 15549676 DOS: 01/21/2009

Flu adjustment appointment. heated out medial side on foot piece to relieve pressure when sitting. Flared pretibial shell. Replaced velcro. Discuss new KAFO. Will submit for authorization and get prescription. will arrange casting afterwards.

Julia

PATIENT CLINICAL INFORMATION

Patient's Name

Practitioner

41998318 NAME: CONSERVA, MAX FABRAZI
DOB: 02/02/1981 M COPAY:
FSC: 743 PLAN: BLUE CROSS/SHLD OF MI 2SU
ATTG MD: UNKNOWN, SCHED 99996
REF MD: BOST, FREDERIC W II MD 77581
AUTH# STOR: OAF LOC: PS9
VISIT: 15277083 DOS: 11/12/2008

copy AFO attached to posterior calf section
allows too much rotation at ankle which
reduces knee control. Max enjoyed pick up of
toe. Switched back to foot section w/ tamarac
will fab new foot piece AFO 3/16" stretched black
copy is solid ankle trimlines. Scheduled next
Friday.

Julie

41998318 NAME: CONSERVA, MAX FABRAZI
DOB: 02/02/1981 M COPAY:
FSC: 743 PLAN: BLUE CROSS/SHLD OF MI 2SU
ATTG MD: UNKNOWN, SCHED 99996
REF MD: NO, REFERRING 99997
AUTH# STOR: OAF LOC: PS9
VISIT: 15287248 DOS: 11/21/2008

11/21/08.
Tried solid ankle AFO piece. instead of
PLS foot piece. Max was happy with it.
Want to do a new KAFO for next year.
will look into that ^{after} ~~over~~ Christmas.

Julie

PATIENT CLINICAL INFORMATION

Patient's Name

Practitioner

41998318 NAME: CONSERVA, MAX FABRAZI
DOB: 02/02/1981 M COPAY: 15
FSC: 743 PLAN: BLUE CROSS/SHLD OF MI 2SU
ATTG MD: UNKNOWN, SCHED 99996
REF MD: NO, REFERRING 99997
AUTH# STOR: OAF LOC: PS9
VISIT: 15059776 DOS: 09/23/2008

Max is breaking medial talar neck joints
size medium in 2 weeks. Too much force
located there. Discussed 2 Danell solutions.
Trying PLS style posterior attachment to
foot piece. Will provide suspension + lateral
control to prevent ankle injuries. Casted over
KATO. Will fit Friday 9/26. A few mods
needed to night splint. Will do ed knee in
2 days.

Julie J

~~9/29/08~~

PATIENT CLINICAL INFORMATION

Patient's Name

Practitioner

41998318 NAME: CONSERVA,MAX FABRAZIC
DOB: 02/02/1981 M COPAY: 15
FSC: 743 PLAN: BLUE CROSS/SHLD OF MI 2SU
ATTG MD: UNKNOWN,SCHED 99996
REF MD: NO,REFERRING 99997
AUTH#
VISIT: 14862816

STOR: OAF LOC: PS9
DOS: 08/06/2008

screws on tamaracks pulled out of. Copoly foot piece. New foot piece was made and fit. Used a longer screw so a nut could be added to prevent screws from backing out. Also gave more gel liners to use as padding. He needs ones that are thicker on one side than the other.

Julif

41998318 NAME: CONSERVA,MAX FABRAZIC
DOB: 02/02/1981 M COPAY:
FSC: 743 PLAN: BLUE CROSS/SHLD OF MI 2SU
ATTG MD: UNKNOWN,SCHED 99996
REF MD: NOT,AVAILABLE 99995
AUTH#
VISIT: 14923380

STOR: OAF LOC: PS9
DOS: 08/22/2008

2 tamaracks (on medial side) ~~APP~~ shredded while playing tennis. Discussed trying w/ longer screws + 2 nuts or going back to old footplate. Max will try at home. If with 2 nuts. At breaks will return for application of old foot plate. Added proximal lateral pad out of Cloud. Interstd in wrap around Night AFO. Submitting for auth. Will see at Mt. Zion on Thursday.

Julif

PATIENT CLINICAL INFORMATION

Patient's Name

Practitioner

41998318 NAME: CONSERVA,MAX FABRAZI
DOB: 02/02/1981 M COPAY:
FSC: 743 PLAN: BLUE CROSS/SHLD OF MI 2SU
ATTG MD: UNKNOWN,SCHED 99996
REF MD: NO,REFERRING 99997
AUTH# STOR: OAF LOC: PS9
VISIT: 14481105 DOS: 05/08/2008

etel liners working great.
Max has added ore at ankle.

lowered medial proximal trimline 3/8"
~~Trimmed foot~~ heated out lateral ankle
to relieve pressure. ON old brace ("swim brace")
cut down foot piece so it was only a strap
and added velcro strap. Flu PRN.

[Signature]
Dorell C.A. Co

1998318 NAME: CONSERVA,MAX FABRAZI
OB: 02/02/1981 M COPAY:
SC: 743 PLAN: BLUE CROSS/SHLD OF MI 2SU
TTG MD: UNKNOWN,SCHED 99996
EF MD: NOT,AVAILABLE 99995
JTH# STOR: OAF LOC: PS9
SIT: 14510077 DOS: 05/12/2008

Carbon fiber brace adjustment.
Heated and flared distal lateral corner on thigh
section, Replaced Strap at ankle and added
teflon strip to ↓ sharpness on edge. Replaced
Stripped screw in ankle joint. Loctited all
screws + provided 3 spares. Flu PRN.

[Signature]
Dorell C.A. Co

PATIENT CLINICAL INFORMATION

Patient's Name

Practitioner

41998318 NAME: CONSERVA,MAX FABRAZI
DOB: 02/02/1981 M COPAY:
FSC: 743 PLAN: BLUE CROSS/SHLD OF MI 2SU
ATTG MD: UNKNOWN,SCHED 99996
REF MD: NO,REFERRING 99997
AUTH# STOR: OAF LOC: PS9
VISIT: 14226781 DOS: 03/05/2008

joint bolts came loose. Reset and re secured uprights on brace

Walt Buffers

41998318 NAME: CONSERVA,MAX FABRAZI
DOB: 02/02/1981 M COPAY:
FSC: 743 PLAN: BLUE CROSS/SHLD OF MI 2SU
ATTG MD: UNKNOWN,SCHED 99996
REF MD: NO,REFERRING 99997
AUTH# STOR: OAF LOC: PS9
VISIT: 14242114 DOS: 03/11/2008

modified footplate with less stiffness than
anaphite worked well in vast majority of cases
would like additional footplate with stiffer

Walt Buffers

plastic for certain activities.

41998318 NAME: CONSERVA,MAX FABRAZI
DOB: 02/02/1981 M COPAY:
FSC: 743 PLAN: BLUE CROSS/SHLD OF MI 2SU
ATTG MD: UNKNOWN,SCHED 99996
REF MD: NO,REFERRING 99997
AUTH# STOR: OAF LOC: PS9
VISIT: 14264325 DOS: 03/17/2008

3/8" copoly footplate. raised ~~lowered~~ medial joint 1/2"
to bring foot flat on ground. DF assist joints didn't
work. He will try it out and return for changes.

Juli J
Doree

41998318 NAME: CONSERVA,MAX FABRAZI
DOB: 02/02/1981 M COPAY: 15
FSC: 743 PLAN: BLUE CROSS/SHLD OF MI 2SU
ATTG MD: UNKNOWN,SCHED 99996
REF MD: BOST,FREDERIC W II MD 77581
AUTH# STOR: OAF LOC: PS9
VISIT: 14458087 DOS: 05/05/2008

Distal medial extension on thigh cuff has too much
pressure. Heated to relieve, removed spenco padding &
placed gel liner section at knee for padding. Used
3mm + 6mm, gel side down. He will experiment and come back
Thursday if needed.

Juli J
Doree

PATIENT CLINICAL INFORMATION

Patient's Name

Practitioner

41998318 NAME: CONSERVA, MAX FABRAZI
DOB: 02/02/1981 M COPAY:
FSC: 743 PLAN: BLUE CROSS/SHLD OF MI 2SU
ATTG MD: UNKNOWN, SCHED 99996
REF MD: NO, REFERRING 99997
AUTH#
VISIT: 14109541
STOR: OAF LOC: PS9
DOS: 02/08/2008


- Orthoplane - Hollywood
818-753-1316

ins.
Jim Dryden, Carmen, Christine

Billing as 2007.


41998318 NAME: CONSERVA, MAX FABRAZI
DOB: 02/02/1981 M COPAY:
FSC: 743 PLAN: BLUE CROSS/SHLD OF MI 2SU
ATTG MD: UNKNOWN, SCHED 99996
REF MD: NO, REFERRING 99997
AUTH#
VISIT: 14126287
STOR: OAF LOC: PS9
DOS: 02/13/2008

Redoing foot shell. Copoly black. Pulling
out over lateral malleoli. Skin intolerant
of pressure thru \bar{c} knee flexion.

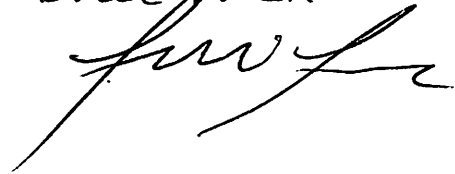

41998318 NAME: CONSERVA, MAX FABRAZI
DOB: 02/02/1981 M COPAY:
FSC: 743 PLAN: BLUE CROSS/SHLD OF MI 2SU
ATTG MD: UNKNOWN, SCHED 99996
REF MD: NOT, AVAILABLE 99995
AUTH#
VISIT: 14164752
STOR: OAF LOC: PS9
DOS: 02/22/2008

Fit \bar{c} black copoly foot shell. Max liked it.
Heated ^{lat.} joint out $\approx 1/8"$. Not sure about pressure
proximal to malleoli. Max will try and let us
know. Also fit protibial shell to distribute
pressure.

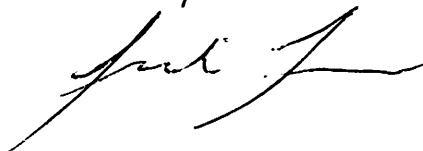

PATIENT CLINICAL INFORMATION

Patient's Name

Practitioner

41998318 NAME: CONSERVA, MAX FABRAZI
DOB: 02/02/1981 M COPAY:
FSC: 1 PLAN:
ATTG MD: UNKNOWN, SCHED 99996
REF MD: NO, REFERRING 99997
AUTH# STOR: OAF LOC: PS9
VISIT: 14052468 DOS: 01/30/2008

black
leather strap 4 1/2" hook
total length

Max is here for fld of double upright. Carbon fiber KAFO - Free knee and ankle. Powder coated bars. Max was comfortable in brace. Holds knee straight. Brace forces him into flexion at about 80° knee flexion. Max wants to try using it for a week.

Copdy KAFO was cracked proximal to ankle joint due to torsion. Walt welded it & top polypro as a temporary hold. Instructed Max not to use new KAFO for snowboarding until we are sure of knee flexion ROM. Will schedule casting for wraparound night splint when he returns from vacation. F/u ~1 week.



2-7-08

max is here for adjustment to KAFO. Added pads where needed. Trimmed medial side 1/2". Added distal calf strap. Added pad to dorsum of foot.

PATIENT CLINICAL INFORMATION

Patient's Name

Practitioner

41998318 NAME: CONSERVA,MAX FABRAZI
DOB: 02/02/1981 M COPAY:
FSC: 1 PLAN:
ATTG MD: UNKNOWN,SCHED 99996
REF MD: NO,REFERRING 99997
AUTH# STOR: OAF LOC: PS9
VISIT: 13846743 DOS: 12/07/2007

✓ osteo. posterior offset creates
problem with axis of motion needs
smaller (thin) steel joints without
posterior offset. Pain joint 1/2". Will put on one joint only
with graftable on other.
Walt Smith MD

41998318 NAME: CONSERVA,MAX FABRAZI
DOB: 02/02/1981 M COPAY:
FSC: 1 PLAN:
ATTG MD: UNKNOWN,SCHED 99996
REF MD: NO,REFERRING 99997
AUTH# STOR: OAF LOC: PS9
VISIT: 13877492 DOS: 12/14/2007

12/14/07
Fit \bar{c} 5/8" x 3/16" stainless steel joint on lateral
side. 1/4" too low. Heat relieved foot plate. Will reevaluate
distal lateral tibial portion after attaching medial joint

41998318 NAME: CONSERVA,MAX FABRAZI
DOB: 02/02/1981 M COPAY:
FSC: 1 PLAN:
ATTG MD: UNKNOWN,SCHED 99996
REF MD: NO,REFERRING 99997
AUTH# STOR: OAF LOC: PS9
VISIT: 13916878 DOS: 12/26/2007

[Signature]

PATIENT CLINICAL INFORMATION

Patient's Name

Practitioner

41998318 NAME: CONSERVA, MAX FABRAZI
DOB: 02/02/1981 M COPAY:
FSC: 1 PLAN:
ATTG MD: UNKNOWN, SCHED 99996
REF MD: NOT AVAILABLE 99995
AUTH# STOR: OAF LOC: PS9
VISIT: 13424300 DOS: 08/24/2007

MAX is a 26yo male here for evaluation and casting for a new (R) KAFO.

Pt seen for casting for new KAFO (R).

Injury ~ 18 years ago resulting in no lateral femoral condyle on lateral tibia. ~ 45° valgus unbraced. Very active, plays tennis and snowboards. Foot pronated but asymptomatic. Will fab check socket and try in 2 wks.

41998318 NAME: CONSERVA, MAX FABRAZI
DOB: 02/02/1981 M COPAY:
FSC: 1 PLAN:
ATTG MD: UNKNOWN, SCHED 99996
REF MD: NO, REFERRING 99997
AUTH# STOR: OAF LOC: PS9
VISIT: 13554131 DOS: 09/26/2007

Julie F
Dental Co. CO.

Initial fitting of (R) KAFO check brace. Fit @ no knee, free ankle joint. Tested WB but no ambulation. Need to pronate foot and remove circumferentially on medial side. Needed changes made. 2nd check socket will be fit after Max returns from Europe.

Julie F
Dental Co. CO.

41998318 NAME: CONSERVA, MAX FABRAZI
DOB: 02/02/1981 M COPAY:
FSC: 1 PLAN:
ATTG MD: UNKNOWN, SCHED 99996
REF MD: NO, REFERRING 99997
AUTH# STOR: OAF LOC: PS9
VISIT: 13570894 DOS: 10/22/2007

2nd check socket fit. Max tolerates corrected pressure and is happy @ corrected position at knee. Will smooth max. Recast then send out for fabrication. Location underevaluated. Will call Max when brace is ready.

Julie F
Dental Co. CO.